

English Version



My Pregnancy Days User Manual


O!fitness **My Pregnancy Days**




Hardware

Front

LCD Screen

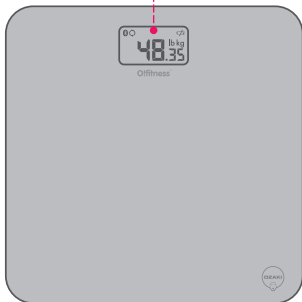


Data sync Low battery
to App

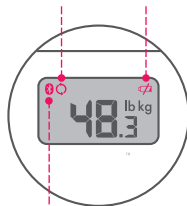

Bluetooth
is connecting

Back

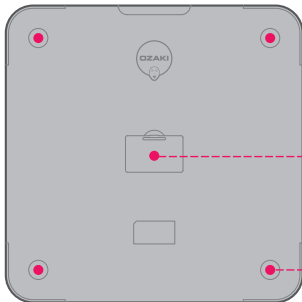


Battery box

Weighing sensors

A: Set up the scale

If you are using the scale for the first time or replacing the batteries, you need to calibrate the scale.

1

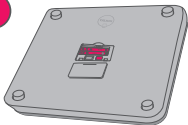


1. Open the battery cover.

2




2. Insert 4 AAA batteries.

3

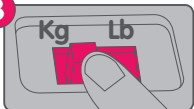

3. Select weight units, Kg / Lb. The button is situated in the battery compartment.



4

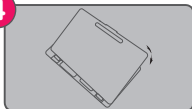

4. Place the cover back.

5. Immediately step on and step off the platform to turn on the scale.

5

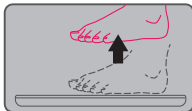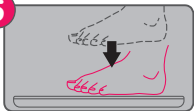

"LCD screen" will show "oooo" then gradually turn into "o", which means the calibration is completed.



B: Pair your scale with your iPhone® /iPad®

1. Download App





or




Scan QR code & download
MyPregnancyDays app

Go to App store & search
MyPregnancyDays & download

2. Turn on the bluetooth

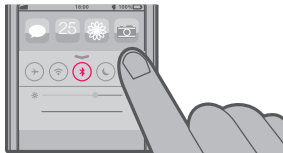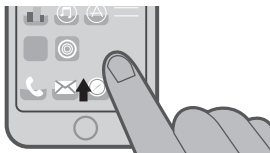

3. Turn on MyPregnancyDays App.

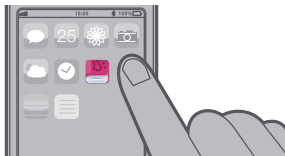

4. Follow the instruction in app to set your personal information.

FAQ & Warming

Something wrong with my scale?

The Scale does not respond

- Check if the batteries are inserted correctly?
- Check the batteries is sufficient or not?

Number keep flashing or the scale cannot measure my weight

- The current weight may be under or beyond the measurement range of 10kg~180kg
- Check if the scale is placed on solid and flat ground.
- Try to stand on the scale stably.
- You should stand on the scale directly. **NOT TAP ON.**

Cannot connect or sync data to my iOS devices

- Check if the bluetooth of your iOS devices is turn on.
- Keep the iOS devices with 2m from the scale
- Make sure Scale is powered on, and bluetooth icon is flashing.

How to clean my scale?

1. Find an old brush and clothes.
2. Brush the surface
3. Wipe the scale with wet clothes.

Anything else I should pay attention to?

- This product is not suitable for pregnant women with twins or more.
- The scale data must be considered purely as reference values and not as a medical opinion.
- Do consult the health care provider if there's any medical issues and follow the advise from your health care provider even it is different from the information My pregnancy days provide.
- Do not use the scale on a wet floor or with wet feet.
- Do not stand on the edges of the scale.



Limited One (1) Year Warranty Information

Warranty Coverage

OZAKI's warranty obligations are limited to the terms set forth below: OZAKI warrants this OZAKI hardware product against defects in design, workmanship and materials under normal use for a period of ONE (1) YEAR from the date of the original user purchase ("Warranty Period"). Accessories accompanying with this product are not covered. This Limited Hardware Warranty may be subject to OZAKI's change at anytime without prior notice.

Repair and replacement

The maximum liability of this OZAKI hardware product shall not exceed the purchase price paid by the original purchaser. The original purchaser shall without undue delay notify OZAKI of any defect which appears according to OZAKI's Policy, failure to which shall mean that the purchaser waives its right to have the defect remedied. A valid form of a bill of sale or receipt from an authorized retailer/distributor with date of the original purchase must be presented to obtain warranty service. If a valid claim is received within the Warranty Period, the sole remedy of the original purchaser and OZAKI's sole and exclusive liability shall be limited to, at OZAKI sole discretion, either repair of the hardware defect using new or refurbished replacement parts, or replacement of the product. Repaired or replacement hardware will be warranted for the remainder of the original Warranty Period or thirty (30) days, whichever is longer. When a product or part is exchanged, the replacement hardware becomes the property of the original purchaser and all hardware or part thereof that is replaced shall become the property of OZAKI.

Exclusions and Limitations

This warranty does not apply: (a) to faulty and improper installation, maintenance, service, repair and/or alteration in any way that is not contemplated in the documentation for the product or canted out with OZAKI' consent in writing, operational adjustments covered in the operating manual for the product or normal maintenance, (b) to cosmetic damages, (c) if the product is modified or tampered with, (d) to damage caused by acts of God, misuse, abuse, negligence, accident, normal wear and tear and deterioration, improper environmental conditions (including, but not limited to, electrical surges, water damage and/or heat/cold/chemical exposure) or lack of responsible care, (e) if the

product has had the model or serial number altered, defaced or removed, (f) to consumables (such as batteries), (g) to products that have been purchased 'as is' and OZAKI, the seller or the liquidator expressly disclaim their warranty obligation pertaining to the product, (h) to any non-OZAKI hardware product or any software (irrespective of packaged or sold with an OZAKI hardware product) and OZAKI products purchased from an unauthorized distributor/reseller, (i) to damage that occurs in shipment or (j) to damages by any other causes not related to defective design, workmanship and/or materials, (k) cracked or broken cabinets, or units damaged by excessive heat.

Warranty Service

The Limited Warranty is applicable in all countries and may be enforced by contacting OZAKI Support worldwide.

For questions and supports, please visit: ozakiverse.com/support

Legal Statement

The information provided along with this product & its APP, either those input by the user and represented by the APP, or related information accessible in the APP or introduced by this product, are popular science information issued from relatively reliable research institution, such as Institute of Medicine and National Research Councils of the National Academies, which must not be deemed or used as a diagnosis for health or a substitute for professional medical assessments or advices. If the User or the pregnant woman that the User records has any questions or any long lasting, pregnancy-related and/or special health issues, please look for consultation and support from professional medical doctors and institutions.

Users of this product do so exclusively at their own risk. This also applies without limitation to any use by the user of data created or provided from the product and its APP, including but not limited to (i) measured data (varied by product functions, e.g. weight, pulse, etc.), or (ii) recommendations for actions (e.g., exercise reminders etc.). The user explicitly acknowledges that any such data or content may contain errors, and Ozaki does not, to the extent permitted by law, assume any responsibility for the correctness of such data. At no circumstance should Ozaki or its associates, dealers or agents, assume the entire liability, cost of repair, correction or any payment or other remedy for the misuse or abuse of this product. Ozaki reserves the right to revise this publication and to make changes from time to time for the contents hereof without obligation or prior notice for the changes.

Specification

Bluetooth connectivity

Battery Requirement: AAA*4 (Batteries not included)

Dimensions: 32.5(W)x 32.5(W)x 2.5 (D)cm

Weight: 2.3kg

Caution & Important Notes

Please read these before use



Do not expose the product to fire or close to high temperature



Stop using if found liquid leakage, overheating or any other unusual conditions



Do not expose the product to water or high humidity



Green power, great recycle using

Customer Service
<http://ozakiverse.com/>



Green for
the Better